



*My* MED

Vendel utca, 11  
1096 Budapest - Hungary  
Tel. +36 70 389 5841

[www.my-med.eu](http://www.my-med.eu)  
[info@my-med.eu](mailto:info@my-med.eu)





# *My* FILLER

NEW LIFE TO THE SKIN



DERMAL FILLERS

# MyFILLER



## PROFESSIONAL, NATURAL, QUALITY

Mr. Roberto Caldeira, the founder of **My Med**, had a vision and dreamt big in 2017. He wanted to create a unique and professional brand of dermal fillers, made of 100% natural ingredients, by taking the concepts of natural beauty and „slow aging” to a higher level. This vision came to life and the success story of **My Med**'s first brand, **My Filler** has been unbroken ever since. The 6 products of the dermal filler brand, have become well-known and loved worldwide, and are present in almost 50 countries on the market of aesthetic fillers by now. A few years after the foundation, in consultation with specialists, we created the **My Filler Revitalize**, another mesotherapy product line of **My Med**. There have been many recognitions and international awards, which have accompanied our journey ever since, for which we are very proud and grateful.

*“Our outer beauty and glow exist and breathe in close unity with our inner, spiritual world, because there is no inner harmony, without outer harmony and vice versa”*

We believe in the power of natural beauty and the importance of increasing self-confidence, and we try to promote this as much as possible with all our products and developments. This is our mission and creed, which we represent on all channels.





## SCIENTIFIC INNOVATION

**My Filler**, the unique family of dermal fillers, includes six different formulas, using two different technologies (**TRIO** and **Dual Link**), all of which are CE certified by an internationally recognized Notified Body in Italy.

In the **Soft**, **Medium** and **Strong** products based on **TRIO Technology**, the specific concentration of hyaluronic acid is combined with different cross-linking rates and molecular sizes, in a special manufacturing environment. As a result, three fillers with largely different properties are obtained, which allows for a wide range of usability.

### High satisfaction on both the client and professional side!

- all desired aesthetic results can be achieved with this composition, nothing is impossible
- high hydrophilic capacity
- easy to shape
- long-lasting effect
- stability after injection



# SOFT




## FINE LINES – GENERAL HYDRATION

The **My Filler Soft** hyaluronic acid filler is recommended for situations where the area to be treated requires a slight increase in volume, the elimination of fine wrinkles, and hydration. It can also be used to brighten the tone of the face, maintain a healthy appearance and delay the formation of wrinkles.

- **Recommended injection depth:** upper dermis
- **Recommended areas:** Periorbital, under-eye circles, perioral, smoker lines, natural-looking – hydrated lips
- **Other characteristics:** extremely low viscosity, so it can also be easily used with various mesotherapy techniques




### Suggested Injection Technique:

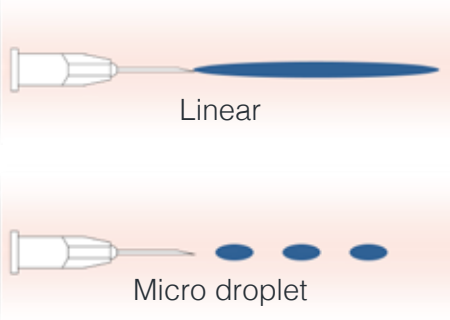

Produced using **TRIO Technology**, it has the smallest HA concentration in the **My Filler** range, just 20mg/ml. It has a very low viscosity making it extremely easy to inject with a 27G or 30G needle.



### Suggested Injection Needle:

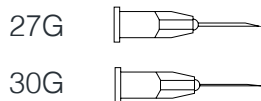


### Suggested Injection Area:



Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF001	20 mg/ml	Monophasic	BDDE	1,5	80

**My Filler Soft pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.

# MEDIUM

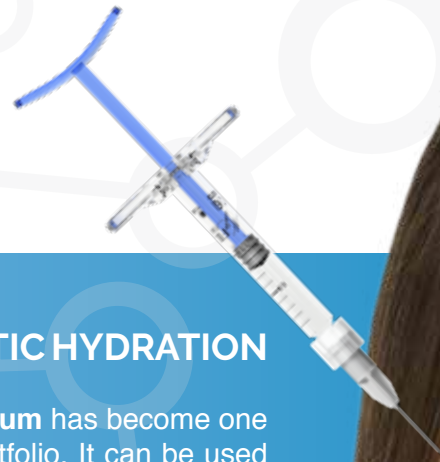

## MEDIUM VOLUME - THE MASTER OF ENERGETIC HYDRATION

Thanks to its universal application properties, **My Filler Medium** has become one of the most popular products from the **My Filler** product portfolio. It can be used perfectly if our goal is to achieve medium volume and energetic hydration in the treated area.




CE 0373

- **Recommended injection depth:** middle part of the dermis
- **Recommended areas:** Nasolabial wrinkles, mentolabial (marionette wrinkles), glabella and lip volume increase
- **Other characteristics:** injecting is easy thanks to the ergonomically designed syringe, the comfortable piston head and the medium viscosity of the material




Produced using **TRIO Technology**, it has 25mg/ml concentration of monophasic HA. It has a medium viscosity suitable for injection with a 27G needle or cannula. Several techniques are suitable with **My Filler Medium** filler such as linear, point or cross-hatching.


### Suggested Injection Area:

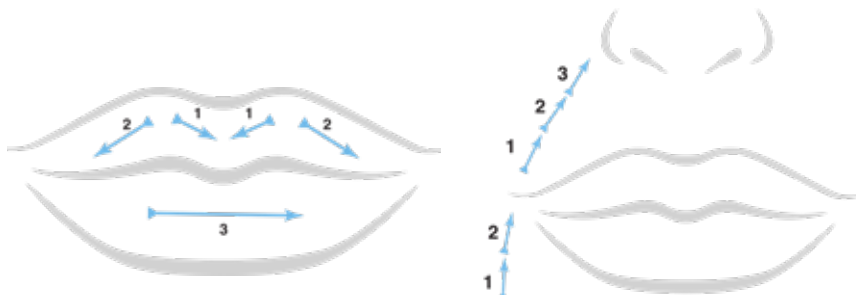

### Suggested Injection Technique:

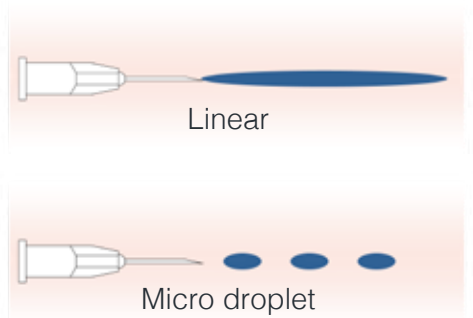

### Suggested Injection Needle:

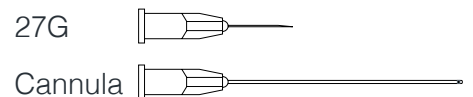

Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF002	25 mg/ml	Monophasic	BDDE	2,5	100

**My Filler Medium pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.

# STRONG

## INNOVATIVE - THE EXPERT IN HIGH VOLUME INCREASE


Due to its high HA concentration, the **My Filler Strong** is unique in the world of hyaluronic acid fillers. Recommended especially for contouring and high volume corrections.

- **Recommended injection depth:** from the deep dermis to the periosteum
- **Recommended areas:** chin, jawline, nose and middle area of the face
- **Other characteristics:** it has an extremely high viscosity and excellent formability properties. After injection, the material remains stable in the area, avoiding any sort of unexpected migration to other areas



### Suggested Injection Technique:

Linear, Micro droplet

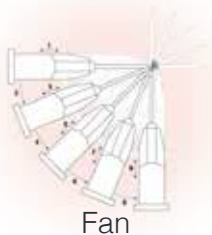


**My Filler Strong** is unique in its type. Produced using **TRIO Technology**, it has the exceptional HA concentration of 30mg/ml. It has a very high viscosity yet has the ability to pass smoothly through needles and cannulas - 25G, 26G and 27G.



### Suggested Injection Needle:

25G - 27G 

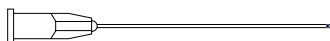Cannula 

### Suggested Injection Area:



Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF003	30 mg/ml	Monophasic	BDDE	3,5	140

**My Filler Strong pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.





## IMMEDIATE AND LONG LASTING RESULTS

It is a well-known fact, that „standard” fillers take longer to achieve the final, desired results. The cross-linked HA needs more time to settle in the desired area and connect to the water molecules. The process takes 5-7 days, depending on the biological characteristics of the patients.

Linear (non-cross-linked, also known as biorevitalizing, skin booster) materials have a very short life span after injection, they are quickly degraded thanks to the body's natural hyaluronidase enzymes. Its effect is that it immediately comes into contact with water molecules, thereby creating a brilliantly hydrated impact, but within 5-7 days linear hyaluronic acid is completely degraded.

The **DUAL Link Technology** was created as a result of years of development and consultations with specialists and patients, which provides a solution to the wishes of most specialists and patients: the achieved result is not only long-lasting but can be seen and felt immediately!

### When the dream of patients and professionals comes true!

The **DUAL Link Technology**, combining the advantages of both linear and cross-linked HA, combines the advantages of both technologies in one unique line of products. The special line using **DUAL Link Technology** provides an immediate effect, thanks to the non-cross-linked technology, as well as a long-lasting result due to the cross-linking HA, ensuring patient satisfaction and saving time for doctors by avoiding future corrections.



# HA Refill




## SKINBOOSTER, ENHANCED TO A NEW LEVEL

The concept of skinboosters is very simple, yet greatly effective. Our experts decided to take the skinbooster to a new level. In addition to the cross-linked hyaluronic acid, non-cross-linked hyaluronic acid was also supplemented with the skinbooster to increase the moisturizing factor. This is how our popular Dual Link technology is connected to the production of this skinbooster, thus the patients can see its immediate effect.


- **The result:** a hydrated, radiant-looking and even-textured complexion, that lasts longer, increased effectiveness
- **Recommended areas:** the entire area of the face
- **Other characteristics:** low viscosity



### Suggested Injection Technique:

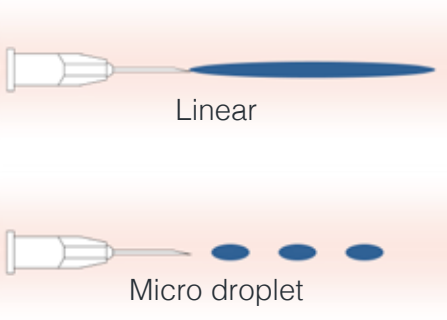


### Suggested Injection Needle:

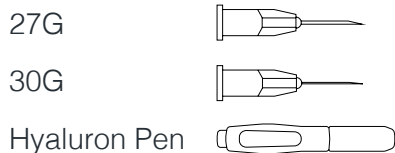

Produced using **DUAL Link Technology**, My Filler HA Refill has 18mg/ml concentration of monophasic cross linked HA and 2mg/ml of monophasic non-cross linked HA. It has a very low viscosity suitable for injection with a 27G or 30G needle or Hyaluron Pen.

It has been proven that monophasic fillers do not migrate within tissues like biphasic fillers. This helps keep excellent results exactly where the client needs them.

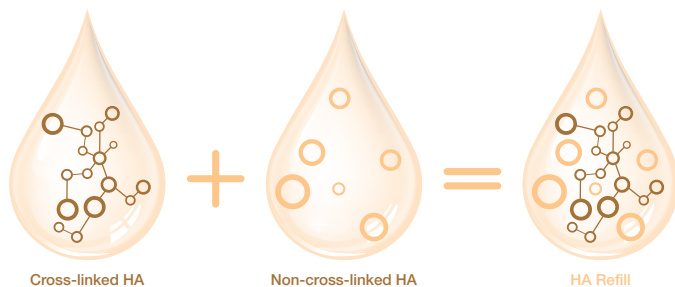

Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF107	18 + 2 mg/ml	Monophasic	BDDE	1,5	80

**My Filler HA Refill pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.

# GLIPS



## LONG-LASTING GLOWING LIPS

**My Filler Glips** is the first product of the **My Filler** range completely dedicated and exclusively made for lip treatments, such as volume and contouring.

This specially composed formula was created with the consideration of the following unique characteristics: biological, color, shape and volume.




- **The result:** immediate and long-lasting shine, full lips
- **Other characteristics:** medium viscosity




CE 0373


Produced using **DUAL Link Technology**, it has 20mg/ml concentration of monophasic cross linked HA and 5mg/ml of monophasic non-cross linked HA. It has a medium viscosity suitable for injection with a 27G needle or cannula. Several injection techniques are suitable with **My Filler Glips** such as linear or point technique.

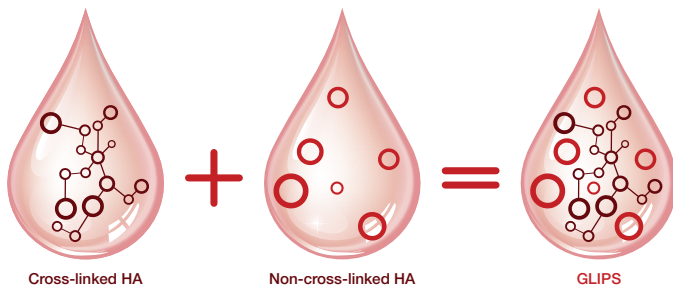

### Suggested Injection Technique:

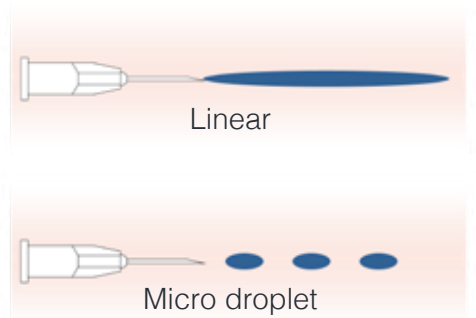

### Suggested Injection Needle:

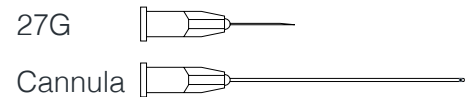

Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF105	20 + 5 mg/ml	Monophasic	BDDE	2	100

**My Filler Glips pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.

# Man Max

## SPECIAL FILLERS FOR MEN


The special formula of **My Filler Man Max** was created taking into account the biological and structural differences between male and female skin. In addition to the strength of facial hair, we can also discover other differences.

- **The result:** a more masculine appearance, a strong jawline and well-defined shape and volume.
- **Recommended injection area:** chin, jaw, nose and midface area




### Suggested Injection Technique:

Linear, Micro droplet

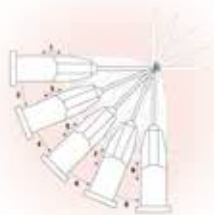


Fan

Produced using **DUAL Link Technology**, **My Filler Man Max** has 25mg/ml concentration of monophasic cross linked HA and 5mg/ml of monophasic non-cross linked HA. It has a high viscosity suitable for injection with a 25G, 26G and 27G needle or cannula. Several injection techniques are suitable with **My Filler Man Max** such as linear, point, fan or cross hatching techniques.

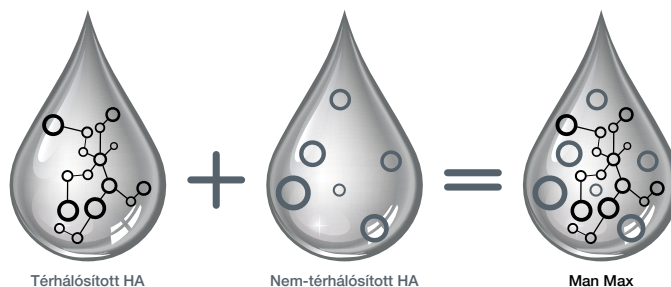

### Suggested Injection Needle:

25G - 27G 

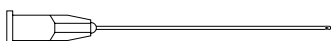Cannula 

Code	HA Concentration	HA Type	Cross-linked agent	Cross-linked rate	Particle Size
MYF004	25 + 5 mg/ml	Monophasic	BDDE	3	120

**My Filler Man Max pack contains:** 1 sterile prefilled syringe with 1 ml, 2 needles 27G 1/2", 2 product labels, 1 multilanguage instruction leaflet, 1 patient card.

## JOIN THE MY FILLER FAMILY

**My Filler** products are available in more than 50 countries around the world. We are constantly working to make our award-winning formulas known in as many places as possible and reach the experts. Our goal is for doctors to enjoy treatments with our high-quality products and to contribute to the satisfaction of their patients by achieving the best results.

Follow us on our social media channels: Facebook, Instagram and YouTube to get more information about the latest news, developments, new products and market trends.

 [www.myfiller.eu](http://www.myfiller.eu)

 [info@my-med.eu](mailto:info@my-med.eu)

 +36 70 389 5841

 [@myfiller](https://www.facebook.com/myfiller)

 [myfillerofficial](https://www.instagram.com/myfillerofficial)

 [MyMedVideos](https://www.youtube.com/MyMedVideos)




## MOBILE APP

**My Filler App PRO** allows you to take “Before” and “After” photographs together with post-treatment recommendations.

An easy to use App which records all your patient information, sets reminders when they are due for more treatment and quickly shares their photos with them.

Practitioners can easily print patient consent forms from the App.

# MyFILLER App Pro

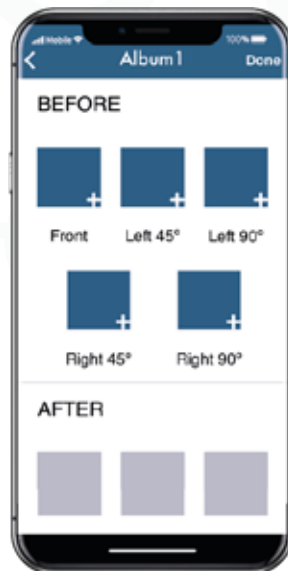

MARKETING INSTRUMENTS

Course Certificate (MYF050)



Clinical Plate (MYF051)



Glips Ice Pack (MYF171)



My Filler Derma Pencil (MYF170)



